

GENDRON[®]

TOTAL
SOLUTIONS
FOR
BARIATRIC
PATIENT
CARE

SAFETY.
COMFORT.
DIGNITY.

The Sentinel II Bariatric Low Bed provides your bariatric patients with all of the above in an easy-to-use design, engineered to allow you to give your patients the best possible care.




Sentinel II Bariatric Low Bed



PATIENT CONTROL PENDANT
Allows patient to operate backrest, knee adjust and nurse call



CAREGIVER CONTROL
Allows caregiver to enable/disable all bed functions from locations on both head-end siderails to prevent patient access; includes electric wheel lock controls

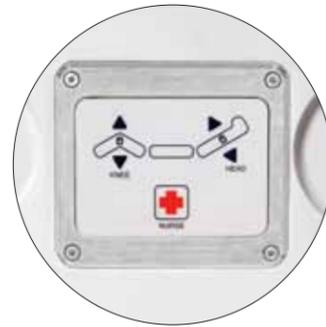

PATIENT CONTROLS
Convenient operation of backrest, knee adjust and nurse call on both head end siderails facing patient




ENCLOSED SIDERAIL MECHANICAL PARTS
Provides additional sanitation and infection control



INTEGRATED BED SCALE SYSTEM
Ensure patient privacy – display can be shielded from view - HIPPA compliant



CAREGIVER CONTROL
Allows caregiver to enable/disable all bed functions from locations on both head-end siderails to prevent patient access

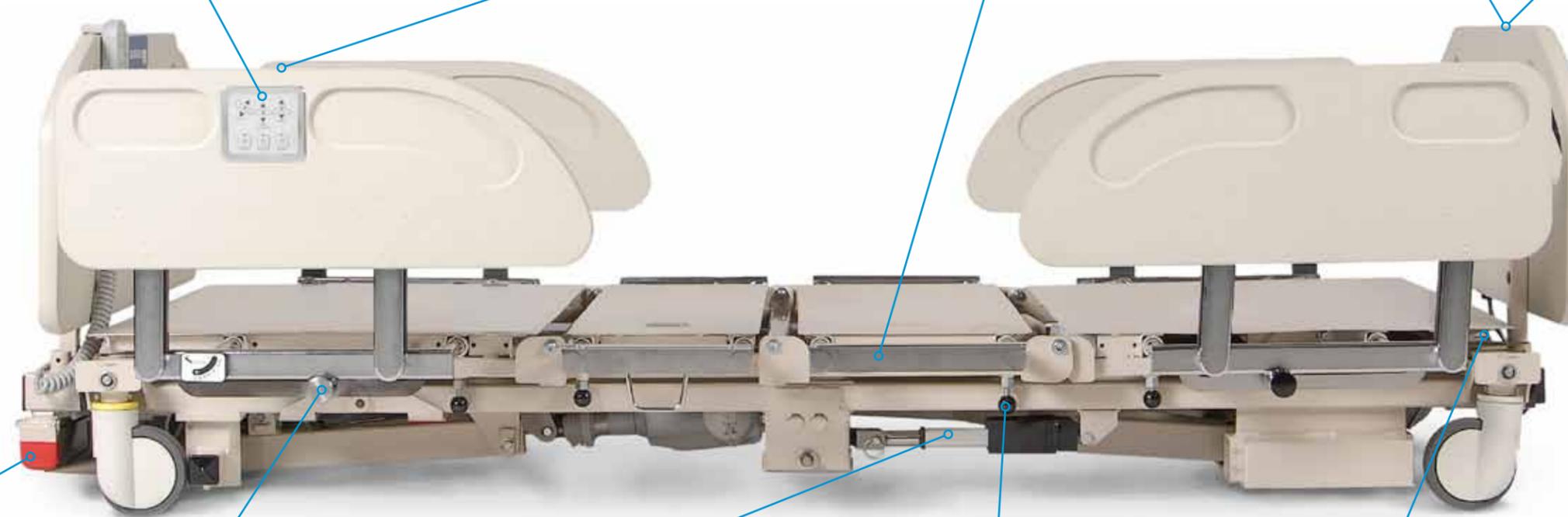

BED POSITIONING/ ARTICULATION
Provide proper clinical positioning from Trendelenburg and Reverse Trendelenburg to knee flex and cardio chair – easy-to-use button controls; one button cardio chair positioning and one touch return to level

OUT OF BED EXIT ALARM
Programmable alarm sounds at the bed and connection to nurse call station (optional)



CPR FOWLER BACKREST RELEASE
Respond faster to patient needs – lowers backrest to flat within six seconds for CPR



EASY-ADJUST SIDE RAILS
Increase patient safety by preventing patient roll over – easy and convenient pull-knob to unlock and lower side rails; 17" high above bed deck



KEYPAD-CONTROLLED 4-WHEEL ELECTRIC CENTRAL LOCK
Ensures all wheels are safely locked with the touch of a button; eliminates trying to locate and access wheel lock pedal



ADJUSTABLE DECK WIDTH
Transport safely and provide comfort – one touch extension from 38" to 48" wide



ADJUSTABLE DECK LENGTH
Ensure the comfort of your tall patients – one-touch extension from 80" to 84" to 88"

SENTINEL II BARIATRIC LOW BED – MODEL 3875 FEATURES / SPECIFICATIONS



FOOT RAISED POSITION



HEAD RAISED, KNEE FLEX POSITION

PATIENT WEIGHT CAPACITY: 850 pounds

FOWLER BACK REST: 0-50 degrees

TRENDELENBURG/REVERSE TRENDELENBURG: 12 degrees

KNEE FLEX/GATCH: 0-20 degrees

DECK LOW POSITION: 9.75"

DECK HIGH POSITION: 28.5"

ADJUSTABLE DECK LENGTH: 80" - 84" - 88"

ADJUSTABLE DECK WIDTH: 38" - 48"

SIDE RAILS (4): 14.45" above deck high position / 2.5" above deck low position

CPR FOWLER BACKREST RELEASE: Head end, left and right

KEYPAD-CONTROLLED 4-WHEEL ELECTRIC CENTRAL LOCK WITH MANUAL EMERGENCY BACKUP:

Wheel lock attendant controls on footboard and both head end siderails

PATIENT HELPER/TRAPEZE: Optional; installs at head end of bed

PATIENT CONTROL PENDANT:

Backrest and knee adjust and nurse call

CAREGIVER CONTROL: Enable/disable all bed functions

BED POSITIONING/ARTICULATION:

High/low, fowler backrest, Trendelenburg/Reverse Trendelenburg, knee flex/leg lift, cardio chair

UNDER BED ILLUMINATION: LED illumination

OUT OF BED EXIT ALARM: Sounds at bed with option to sound at nurse call station using DB37 cable or 1/4" phone connector

INTEGRAL SCALE WEIGHING SYSTEM:

- Scale display rotates – HIPAA-compliant
- Powered by bed electrical system and battery backup – does not require batteries and will function when disconnected from A/C power
- Weight change mode
- Dual weight display
- Freeze mode
- Backlit display
- Accurate to within 0.2% or 2 pounds, whichever is greater
- Scale displays accurate weight when bed is in ANY articulated position and ANY caster orientation
- 120V power outlet at foot end of bed

WHY BUY GENDRON:

- 26 years of bariatric experience
- Made in the USA
- Complete bariatric product line

GENDRON 

800.537.2521 • sales@gendroninc.com

www.gendroninc.com